

NASAL ENDOTRACHEAL TUBE TAPING



Step 1. Clean the peri-nasal area and dry. Apply benzoin tincture to all areas to be covered by elastoplast



Step 2. Place first tape at the edge of the nare on **the opposite side** to the ET tube



Step 3. Bring the bottom leg straight across under the nose, avoiding the lip



Step 4. Place the top leg over the nose, down and around the ETT with 2 tight anchoring turns




Step 5. After two anchoring turns, wind the tape up the ETT by 1-3cm



Step 6. Place 2nd tape at corner of the nare the ETT passes through



Step 7. Place the upper leg up over the nose and adhere over the first tape

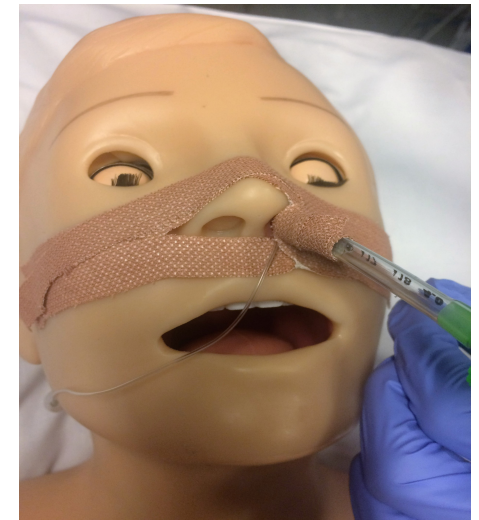

Step 8. Bring the lower leg up around the ETT. After 2 anchoring turns, wind up the ETT and leave small tab for removal